


NEXT-GENERATION AIRWAY TOOLS-DESIGNED WITH THE VIDEO LARYNGOSCOPE IN MIND

Airway tools designed to take full advantage of the superior view provided by video laryngoscopes. The Through-The-Cords Articulating Introducer allows precise tip control. Designed to make the most difficult airways easy. Once in the trachea, an intuitive depth control system allows precise tip depth control throughout the intubation.

Total Control Introducer™ (TCI Si)

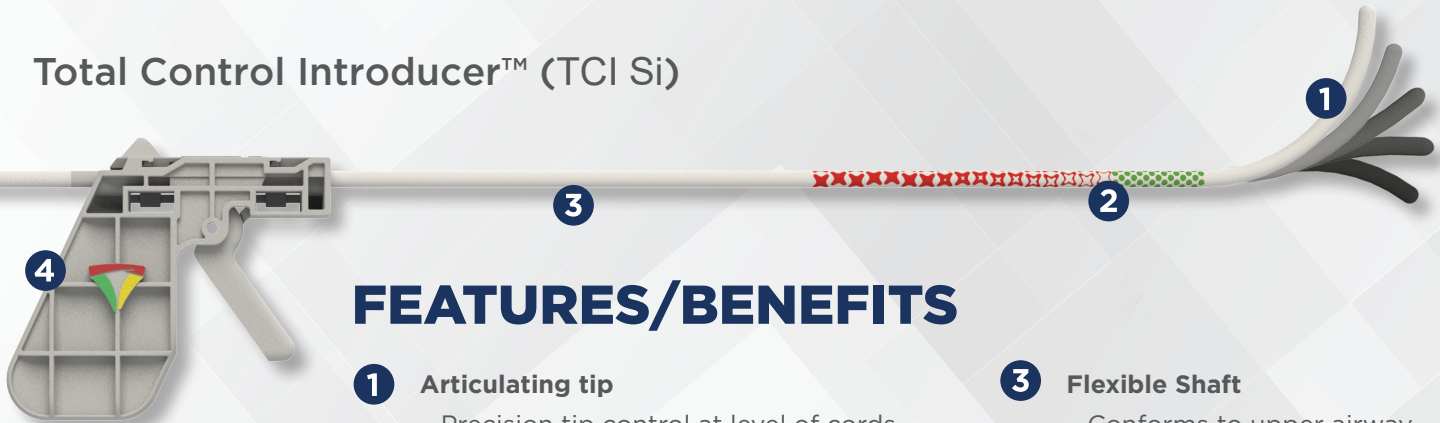

FEATURES/BENEFITS


- | | |
|--|---|
| <p>1 Articulating tip</p> <ul style="list-style-type: none"> - Precision tip control at level of cords - Easy access into the trachea | <p>3 Flexible Shaft</p> <ul style="list-style-type: none"> - Conforms to upper airway - Suitable for nasal intubations |
| <p>2 Intuitive Depth Control System</p> <ul style="list-style-type: none"> - Designed to prevent injury in lower airways - Keeps tip safely mid-trachea | <p>4 Detachable Handle</p> <ul style="list-style-type: none"> - No need to preload tube before use - Ready for immediate use |

ADDITIONAL TOOLS

Intuitive Depth Bougie™ (IDB)



Intuitive Depth Exchange™ Catheter (IDE)

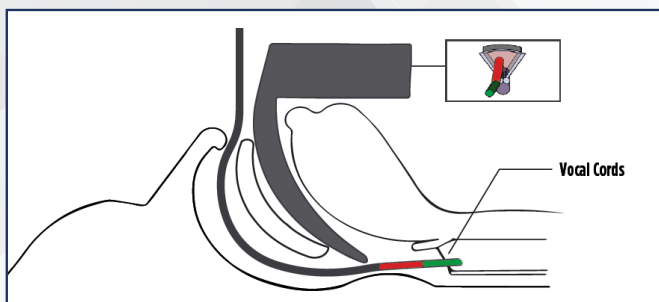

If you can see the cords, you can intubate the cords. Green at the cords, tip mid-trachea.
Simple as that.

TOOLS DESIGNED TO REDUCE FIRST PASS FAILURE AND POOR DEPTH CONTROL



Intubation: The Fastest Way to Save a Life (or End It)

- We can only live for **3-5 minutes without oxygen**
- **40M intubations** per year in the USA (4M ED and ICU intubations)

FIRST ATTEMPT FAILURE, POOR DEPTH CONTROL AND THE COST OF CURRENT COMPLICATIONS

600K → \$20B → 120K

Major Complications

Complications Cost

Deaths

Product Name	Description	Quantity	Product Code
Total Control Introducer™ (TCI Si)	Articulating Introducer with detachable handle, intuitive depth system, 70 CM, 15 Fr, latex free, single use, sterile, fits 6.0 and larger ID tubes	5/Box 20/Case	000946
Intuitive Depth Bougie™ (IDB), Coude Tip	Bougie with intuitive depth system, 70 CM, 15 Fr, latex free, single use, sterile, fits 6.0 and larger ID tubes	10/Box 100/Case	000118
Intuitive Depth Exchange™ Catheter (IDE)	Catheter with intuitive depth system, 70 CM, 15 Fr, latex free, single use, sterile, fits 6.0 and larger ID tubes	10/Box 100/Case	000294

CONTACT INFORMATION

Phone: 1 (801) 917-3851 | www.TTCmed.com